

Hair clippers, nail clippers, suction, lights, tables:
All the bits and pieces that are part of everyday
veterinary practice but don't quite fit anywhere else.

We've selected the items in this section over many
years. If you can't see what you need, please call us. We
may have hidden it somewhere else. If not, we may be
able to source for you.

We offer a wide range of training aids to support you.
From books, posters to step by step guides and
anatomical models.

Custom Furniture

Vet Sink

Major Vetsink features include -



- Fully welded internal joints precluding matter ingress and the retention of infectious organisms.
- All round accessibility in standing and seated positions.
- Sloping sink base draining all fluids automatically away from the subject.
- Integral mixer taps and umbilical shower unit.
- Easily cleaned and maintained, thus popular with veterinary nurses.
- Radiused corners.
- Can be combined with the latest air driven dental units and mobile seating.
- Four button hooks are located on the sides for hanging the taps during cleaning.
- 660mm width.

VETSINK

VETSINK12 Vet Sink 1200mm Tub Table complete*

VETSINK15 Vet Sink 1500mm Tub Table complete*

VETSINK20 Vet Sink 2000mm Tub Table complete*

VETSINKTOP Solid Top converts above into Exam Table/ Op Table*

*Carriage is extra for all items marked **

Tables

Stainless Steel Furniture can be fabricated to your design in our Sheffield Workshop. **Phone or fax a sketch for a quote.**

Folding Wall Mount Examination Table



Very rigid and will not tip over. Surface area can therefore be smaller than a free standing table.

Table flips flat against the wall for cleaning and examination of large dogs.

Standard Table is 40" x 20", but all sizes available on request.

FOLDING WALL MOUNT EXAMINATION TABLE

028026 Folding Table 40" x 20" *


*Carriage is extra for all items marked **

Scrub Dispenser Stainless Steel Demountable

Scrub Dispenser - Demountable



025134 Dispenser

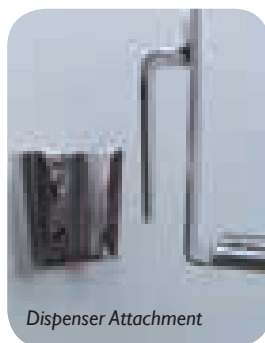

Dispenser Attachment

- Stainless construction.
- Unit dismountable from wall bracket (included) for cleaning/autoclaving.
- Operates standard scrub bottles.

SCRUB DISPENSER - DEMOUNTABLE

025134 Scrub Dispenser Stainless Steel

Mayo Instrument Table



- Stainless Steel throughout.
- Tray 19" (490mm) X 13" (330mm).
- 4 castors.
- Height adjustable 33" (840mm) to 47" (1200mm).
- Lifetime guarantee.

MAYO INSTRUMENT TABLE

028027 Stainless Mayo Instrument Table 19" x 13"*

*Carriage is extra for all items marked **

Scales

Small Platform Scales



New Stainless steel compact Scales, for lab or small animals

SMALL PLATFORM SCALES

025131 Small Platform Scale

Medium Platform Scales



For patients weighing between 1Kg and 20Kg this Medium Table Scale is ideal. Mains and battery power 520mm x 200mm.

MEDIUM PLATFORM SCALES

025135 Medium Platform Scale

Henry Schein Walk On Scales



Developed in Germany exclusively for Henry Schein Animal Health, with vets in mind.

These Scales are robust, durable and easy to use, with a unique easy to clean groove.

Features:

Clear, bright LCD display made of high quality impact proved plastic.

Mains or battery powered – mains adaptor included

Platform size :900x600mm

Low profile platform at only 40mm high

200kg maximum weight

HENRY SCHEIN WALK ON SCALES

025136 HS Walk On Scales

Microscope

BioBlue Microscope



The BioBlue Microscope has been very popular across Europe with veterinarians, and Veterinary Instrumentation is pleased to supply it in the U.K. This high quality Binocular Microscope offers high quality at the right price.

Technical information:

- High quality DIN WF 10x/ 18 Eyepieces
- Standard Objectives 4x, 10x, S40x and S100x (oil) giving total magnification of 40x, 100x 400x and 1000x
- LED 1 watt adjustable illumination which generates no heat. NB Always start with the light on low and increase to level required. This light is much brighter than conventional bulbs!
- Sensitive Focus Control. Coaxial coarse and fine adjustment knobs on both sides with 0.0020mm graduations.
- Stable ergonomic design
- Supplied with Dust Cover and free Aluminium Storage Case.

BIOBLUE MICROSCOPE


BB4260 BioBlue Microscope

Lighting

Good lighting is an essential pre-requisite for any examination or procedure. Regardless of the quality and range of your surgical instruments if you cannot adequately see what you are doing, your procedure will be compromised.

We have selected a range of quality lights to meet your examination and procedure requirements.

Magnifier

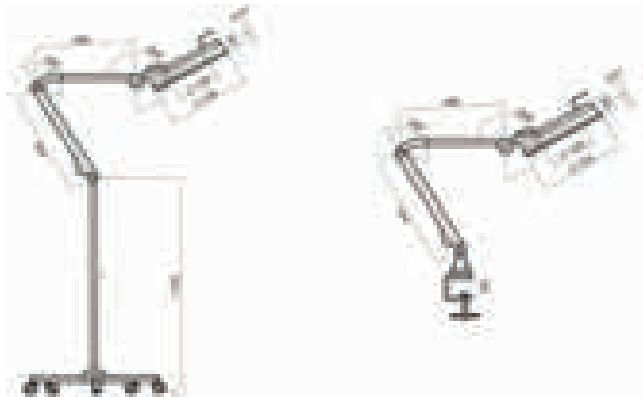

A superior Magnifying Light for detailed distortion free examination.

- 2.25x magnification (5 dioptr).
- 127mm diameter.
- 3 adjustable joints.
- Fluorescent tube 22w.
- 230-240v 50Hz AC.
- 10,500 lux light intensity.
- Options:
 - C-Clamp 6.8kg
 - Mobile 10.0kg

MAGNIFYING PROCEDURE LIGHT

- RLMCLT** Mobile *
- RLMCL** Mobile with 'C' Clamp *

Carriage is extra for all items marked *

HF Series Multipurpose Examination Lights



12929

12930

12454



12455 (Clamp)

12456 (Wall Mount)

This very compact single light source is an ideal procedure light. The head swivels around to whatever position you require.

- Vented shade.
- stay cool finger tip controller.
- 50w halogen light source.
- 230-240v 50Hz AC.
- 31,000 lux light intensity.
- Options:
 - Articulated Wall Mount
 - Articulated Ceiling Mount
 - Flexible Goose Neck Mobile
 - Flexible Goose neck 'C' Clamp



HF SERIES PROCEDURE LIGHT

- 12456** HF Flexible Goose Neck Wall Mount 50w *
- 12929** HF 3-Part Ceiling Mount 50w *
- 12454** HF Flexible Arm Mobile Base 50w *
- 12455** HF Flexible Arm 'C' Clamp 50w *
- 12930** HF Articulated Wall Mount Procedure Light *

Carriage is extra for all items marked *

Due to various ceiling configurations - ceiling lights are not supplied with mounts - please contact info@vetinst.com for further details.

TriStar

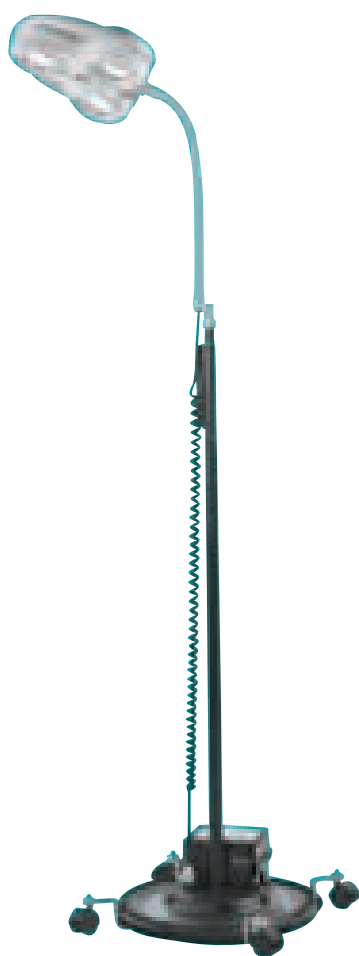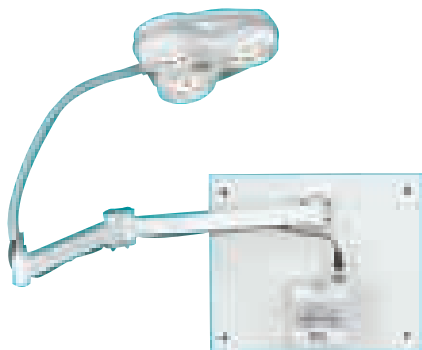

A compact combination of three light sources offers the following features:



- Light pattern is 125mm at 600mm distance.
- 3 x 20w long life halogen bulbs (2000hr).
- 230-240v 50Hz AC.
- 42,000 lux light intensity.
- Light colour 4,100 Kelvin.
- Wall mount option offers 1400mm reach. 8.2kg.
- Mobile unit raises from 685mm to 2030mm on a non tip base. 11.4kg.

TRISTAR LIGHTING

TS2015W	Wall Mount TriStar *
TS2015M	Mobile TriStar *
I8464	Bulb 20w

Carriage is extra for all items marked *

Celestial Star

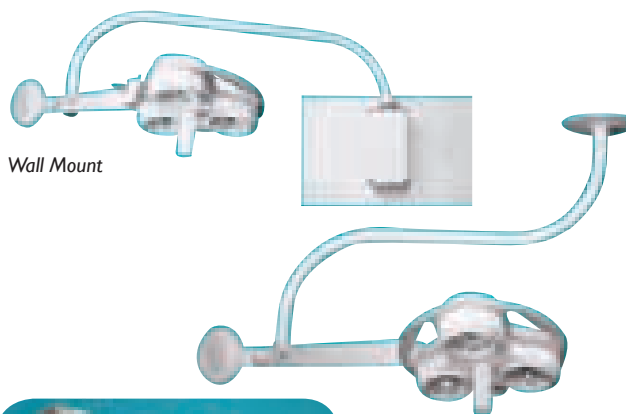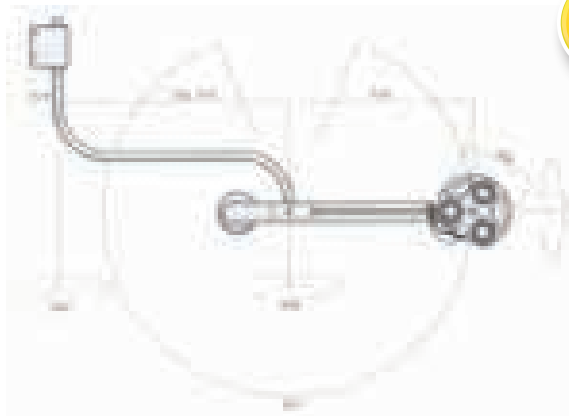

Wall Mount

Ceiling Mount



Mobile CS Unit

Compact and counterbalanced, the three light source Celestial light is very easy to move and adjust.

- Drift free balance arm.
- 3 x 35w long life halogen bulbs (4000hr).
- 230-240v 50Hz AC.
- 64,500 lux light intensity.
- Light colour 4,100 Kelvin.
- Options:

- Ceiling 23.6kg
- Dual Ceiling 47.2kg
- Wall 17.3kg
- Mobile 24.3kg

CELESTIAL STAR LIGHTING

CS2050C	Ceiling *
CS2050W	Wall Mount *
CS2050M	Mobile *
I8544	Autoclaveable handle
I3078	Bulb
I8676	Corrective Filter

Available to UK only. Carriage is extra for all items marked *

Suction Units

HospiVac 350



The HospiVac 350 has a greater maximum vacuum (-0.9Bar) greater flow rate (60 litres/ min) and a greater capacity (2 x 2 litre bottles) than our standard Suction Unit. It is suitable for hospital situations where greater volumes of fluid might expect to be aspirated. The double bottle facility further increases capacity. Size 600 x 460 x 420mm.

A feature found on the HospiVac 350 which although desirable is not found on base level suction units is the foot control. Without this feature the unit is essentially on or off all of the time.

HOSPIVAC 350

185H350	Hospivac 350 Complete with Foot Control *
185H210351	Spare Bottle
185H210353	Spare Bottle with Lid
185HSP033	Spare Tubes
185HSP0121	Spare Filters
185HSP0032	Spare Tubes and Filters

* Carriage is extra for all items marked.

Askir Suction Machine



This simple Suction Unit is sensibly priced to make it suitable for general practice as well as more specialised centres.

The Suction Tube will fit all Standard Tips.

The Bottle capacity is 1 litre which is adequate for virtually all small animal procedures.

The in-line filtration system minimises biohazards.

Most veterinary dentists now regard suction as being essential.

Suction is controllable and there is a vacuum indicator.

CE marked for peace of mind and quality.

The Suction Machine is ready for the attachment of tubes and tips.

ASKIR SUCTION DEVICE

185000	ASKIR 30 Suction Device
185003	Filter Only
185004	Spare Suction Bottle & Lid
185005	Spare Lid from Suction Bottle

Suction Tips

We stock the four most useful styles of tip. Each Tip comes complete with a length of Autoclavable Silicone Tube to connect to the suction unit.

Yankauer



General purpose suction tip for removal of large volumes e.g. from abdomen.

YANKAUER SUCTION TIP

1854327	Yankauer Suction Tip
----------------	----------------------

Poole



As Yankauer but more holes for greater volume uptake. With perforated outer sheath removed the Poole functions as a small bore straight suction tip.

POOLE SUCTION TIP

1854330	Poole Suction Tip
----------------	-------------------

Adson 4mm




Focal suction tip removes fluid from specific locations e.g.in spinal procedures.Very controllable at the handpiece by single finger over air bypass hole.

ADSON SUCTION TIP

I85212I Adson Suction Tip 4.0mm

Frazier 2.6mm

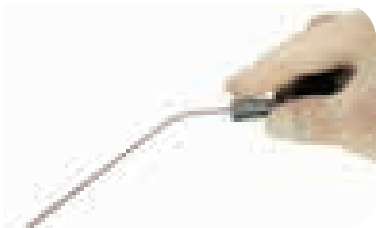

As the Adson but finer still.

FRAZIER SUCTION TIP

I854328 Frazier Suction Tip 2.6mm

Disposable Suction Tips and Autoclavable Tip Holder



The Autoclavable Suction Tip Holder fits onto Standard Silicone Suction Tubing. There are a wide range of tips which fit the standard 11mm tapered opening. We offer a choice of two, narrow and wide. The autoclaveable nature means any standard suction unit may be used for both dentistry and sterile surgery.



DISPOSABLE SUCTION TIP


- I85DISP** Disposable Tips Narrow Mouth (Pack of 3)
- I85DISPW** Disposable Tips Wide Mouth (Pack of 100)
- I85DISPH** Disposable Tip Handpiece (Autoclavable) to fit onto Standard Suction Tubing

Suction Tubing

SUCTION TUBING

I85001 Autoclavable Suction Tubing per Metre

Suction Tubing Connector

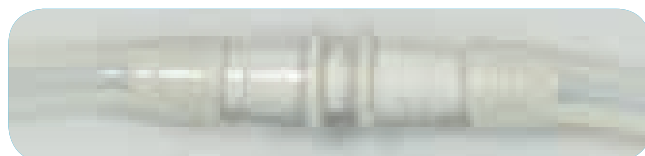

SUCTION TUBING CONNECTOR

I85006 Suction Tubing Connector

Cautery Units

Portable Cordless Cautery Kit



This handy Set comprises a Handpiece using two AA batteries (included) and four interchangeable and replaceable Cautery Tips

- 1/2"(12mm) Fine general purpose electrode for haemostasis
- 1/2"(12mm) Loop electrode for small skin tumours
- 4"(130mm) Fine electrode for deep areas
- 4"(130mm) Loop electrode

Tip temperature 2200°F/1200°C


PORTABLE CORDLESS CAUTERY KIT

- 02313** Portable Cordless Cautery Kit
- 02313W** Replacement Cautery Tip 1/2" Fine
- 02313X** Replacement Cautery Tip 1/2" Loop
- 02313Y** Replacement Cautery Tip 4" Fine
- 02313Z** Replacement Cautery Tip 4" Loop

Thermocautery Unit



This Mains (240v) Thermocautery Device is very robust and reliable with a wide selection of tips. As the trigger is pressed the area to be cauterised is illuminated by an integral light. Suitable for a wide range of minor procedures where thermocautery is appropriate.

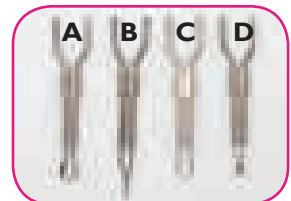

THERMOCAUTERY UNIT

- SP33** Thermocautery Unit 240v
- SP36** Thermocautery Tip A Knife
- SP34** Thermocautery Tip B Pin
- SP35** Thermocautery Tip C Loop
- SP37** Thermocautery Tip D Ball
- SP43** Spare Bulb

Electrosurgery

GIMA BiPolar Electrosurgery Unit 120 watt



HF 122 is a complete innovation in the electrosurgery field. It is a high frequency-electronic scalpel suitable for microcoagulation, small and medium MonoPolar surgery and BiPolar coagulation. The following MonoPolar modes of operation are available via the control panel: pure cutting (CUT), coagulated cutting (COAG), MonoPolar coagulation and microcoagulation (BLEND). Additionally with BiPolar forceps, bipolar coagulation and microcoagulation are available.

Safety Alarm

Open Circuit (OC) Alarm warns the user that connection into the neutral plate has not been made or that the cable is disconnected. The alarm is signalled with an intermittent beep together with a flashing red light on the front panel. In this condition HF emission for cut and coagulation operation is switched off. Meets all EC standards for this type of equipment

Overheat Alarm

Should the final module in the HF generator overheat, the unit incorporates an internal circuit protection which cuts off operation. In this condition the red OVH warning light illuminates and the user must wait until the generator cools down, about 20 sec, before recommencing. The module will not go into protection mode as long as the instruction for operation as described in the manual are followed.

Kit includes:

- Two button MonoPolar handle (Autoclavable)
- Two button MonoPolar handle (disposable)
- Flat neutral plate
- BiPolar cable and forceps
- Foot control
- 6 Short Mono Polar Electrodes



6 x MonoPolar Tips included in the Kit.
Spare Tips available.



BiPolar Forceps are included in the set.

Other style forceps are available.

EU Forceps have a square fitting into the cable.

GIMA BIPOLAR ELECTROCAUTERY 120 WATT UNIT

- HF122** 120 watt MonoPolar/ BiPolar Electrosurgery Unit
- HF125** Spare Autoclavable MonoPolar Handpiece for HF122
- HF30547** Spare MonoPolar Handle Screw Grip for HF122
- HF123** Spare Standard Bipolar Forceps for HF122
- HF124** Spare Bipolar Cable for HF122 (A Type)
- HF30568** Spare Foot Pedal for HF122
- HF127** Waterproof Foot Pedal for HF122

GIMA BiPolar Electrosurgery Unit 160 watt



This Unit is a higher-powered version of our popular HF 122 Unit.

The additional power availability comes into its own when working with mixed tissue densities such as in orthopaedics, and in more moist situations such as mucous membranes.

Modes of operation are pure cut, coagulated cut, forced coagulation, soft coagulation, and with bipolar forceps bipolar coagulation.

HF160 has the same set of alarms as the 120, including Open Circuit which warns of a disconnected cable or poor patient plate contact and interrupts the HF emission for cut and coagulation.

All the standard accessories shown for the HF120 fit this unit. We can also supply an adaptor which allows standard 2 x 4mm banana plug (Valley Lab type) bipolar cables to be used with HF 160.

Technical specifications:

Teil	Description	HF 120	HF 160
-	Electrosurgical unit cable COAG	0000	0000
-	Electrosurgical unit cable CUT	0000000000	0000000000
1 (2%)	Minimum powerable power	0	0
-	Level stop	1	1
-	Digital level display	■	■
1 (2%)	Minimum output power CUT (W)	10 - 200	10 - 200
1 (2%)	Minimum output power BLEND (W)	10 - 200	10 - 200
1 (2%)	Minimum output power COAG FORCED (W)	10 - 160	10 - 160
1 (2%)	Minimum output power COAG SOFT (W)	10 - 160	10 - 160
1 (2%)	Minimum output power BIPOAR (W)	10 - 160	10 - 160
1 (2%)	Modulation factor CUT	Fast 100%	Fast 100%
1 (2%)	Modulation factor BLEND	Fast 100%	Fast 100%
1 (2%)	Modulation factor COAG FORCED	Fast 100%	Fast 100%
1 (2%)	Modulation factor COAG SOFT	Fast 100%	Fast 100%
1 (2%)	Modulation factor BIPOAR	Fast 100%	Fast 100%
1 (2%)	Coag Factor CUT	1:1	1:1
1 (2%)	Coag Factor BLEND	1:1	1:1
1 (2%)	Coag Factor COAG FORCED	1:1	1:1
1 (2%)	Coag Factor COAG SOFT	1:1	1:1
1 (2%)	Coag Factor BIPOAR	1:1	1:1
1 (2%)	Working frequency	500 kHz	500 kHz

GIMA BIPOLAR ELECTROCAUTERY 160 WATT UNIT

- HF160** 160 watt MonoPolar/ BiPolar Electrosurgery Unit
- HF30642** BiPolar Cable Adaptor for Gima Unit HF160 to 2 x 4mm Banana Plus

BiPolar Dissecting Forceps



BiPolar Cable Generator Terminations



A GIMA HF122
B Storz
D 4mm Banana Plugs
 Valley Lab
 Eschman
 Arthrex OPES


These are the most common types of termination for the generator or machine end of the cable. We can supply cables to fit all machines. If you do not see the right connector here please check the brand and model number of your electrosurgery unit and call or e-mail us. We will find the right cable for your machine.

Our Bipolar Cables and forceps all have an EU fitting square end.

BIPOLAR ACCESSORIES

- HF190-100** Bipolar Forceps Adson 1mm Tip
 - HF120-100** Bipolar Forceps Straight 1mm Tip 150mm Long
 - HF130-100** Bipolar Forceps Straight 1mm Tip 180mm Long
 - HF122-100** Bipolar Forceps Straight/ Angled 1mm Tip 150mm Long
 - HF132-100** Bipolar Forceps Straight/ Angled 1mm Tip 180mm Long
 - HF124** Spare Bipolar Cable A Type EU Fitting Forcep Connector
 - HFU-30** Spare Bipolar Cable B Type EU Fitting Forcep Connector
 - HFS-30** Spare Bipolar Cable D Type EU Fitting Forcep Connector
- Other Bipolar Cable Types available. Phone for details

GIMA MonoPolar Electrosurgery Unit

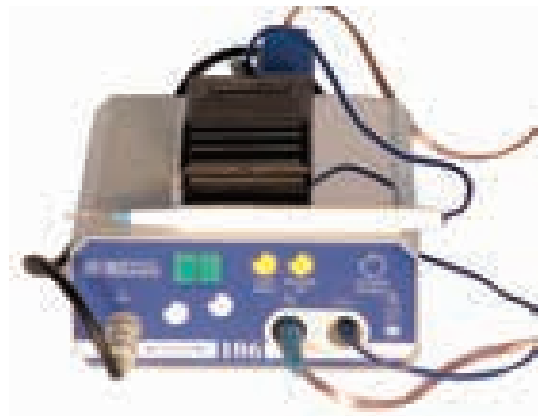

This HF 106 Electrosurgery Unit is particularly suitable for small procedures. The 50 watt unit offers precision of cut with immediate coagulation. There is a choice of 10 electrodes. The three modes of operation, CUT, COAG and BLEND, can be easily selected through a switch on the front panel. Meets all EC standards for this type of equipment. The package includes HF 106 with all accessories including electrode handle, flat neutral plate, foot switch and a kit of 6 electrodes 5 cm long.



6 x MonoPolar Tips included in the Kit.
Spare Tips available.

MONOPOLAR ELECTROSURGERY UNIT

- HF106** 50watt MonoPolar Electrosurgery Unit
- HF107** Spare Autoclavable MonoPolar Handle for HF106

MonoPolar Spares and Accessories



MONOPOLAR ACCESSORIES

- HF1** Electrode
- HF2** Electrode
- HF3** Electrode
- HF4** Electrode
- HF5** Electrode
- HF6** Electrode
- HF7** Electrode
- HF8** Electrode
- HF9** Electrode
- HF10** Electrode
- HF1-10** Electrodes Mixed Set of 10
- HF30537** Scratch Pad (supplied sterile) Pack of 10
- HF30563** Spare Cable for Gima Plate
- HF126** Spare Metal Plate for HF122 & HF106

Hyfrecator® 2000 Electrosurgical System




The ConMed Hyfrecator® 2000 is a multi-function electrosurgical generator with both monopolar and bipolar modes, self-diagnostics, and recall of stored power settings.

Simple and safe to use, the Hyfrecator® has dual microprocessor control. One processor controls the output and functionality of the machine, while the second monitors the system and terminates activity if a safety hazard occurs for both patient and operator safety.

The Hyfrecator® has 2 modes of output. Monopolar has 2 settings, Hi and Lo, and a Bipolar option is also available. The construction of the unit only allows one option to be active at any time, mechanically locking out the connectors for the unused options for added safety.

The independent outputs produce two different therapeutic effects, giving better control of the desired clinical outcome. Lo gives fine control in 1/10 watt increments, while Hi gives longer, broader arcs of current.

The monopolar Lo setting has the added advantage of not requiring the use of a patient plate, which speeds up patient setup considerably as there is no need to clip the plate site as usual. The patient only needs to be in good contact with the table. This setting is suitable for light desiccation and fulguration procedures requiring low intensity and increased precision.

Hi setting can be used with or without the patient plate, but for coagulation whichever power level is used a patient plate must be used, as this improves the coagulation effect by providing a better path for the current flow.

Good wide contact with the patient plate (also called the dispersive plate) is essential to avoid burns, as with all electro-surgery units. The contact area should be clipped, and should be a clean, dry, well-vascularised muscle mass. Bony prominences should not be in contact with the plate.

Monopolar handles are available with fingerswitch controls which allow current settings to be changed from the handpiece, or with no switches which are controlled by an optional footswitch.

Unit Specifications

Depth: 102mm
 Width: 222mm
 Height: 190mm
 Weight: 2.7kg

Output power:

Low: 0-20 watts (1/10 watt increments up to 10 watts)

High: 0-35 watts

Bipolar: 0-35 watts



Hyfrecator® 2000 should be wall-hung or mounted on this convenient optional rolling stand.

What do the terms mean?

Monopolar A single electrode in a pencil handpiece passes current through the patient to a plate electrode which returns the current to the unit. If no plate is used, as with the Hyfrecator®, current returns to ground via good contact with the table.

Bipolar Insulated Bipolar forceps are used to grasp tissue and current is applied using a footswitch – current passes between the tips of the forceps through the grasped tissue.

Desiccation 'drying' of tissue by evaporating the cellular fluids, blanching the treated area. Effect increased by time/power and size of electrode.

Fulguration as desiccation but the electrode is held away from the tissue, resulting in sparking to the surface. Allows a broader sweep of tissue to be treated quickly.

Coagulation causes blood and tissue to clot. Monopolar coagulation with a plate causes heavy dehydration at the electrode site. Bipolar coagulation is only between the tips of the forceps so there is much less collateral tissue damage.

Setup Variants and Accessories

Hyfrecator® 2000 is available as a Monopolar setup with Fingerswitch Pencil Handpiece, Monopolar with Footswitch or complete Mono and Bipolar set-up.

Monopolar Pencils

These are autoclavable, re-usable pencils handpieces.

The Fingerswitch Pencil Hand Control has with power adjustment buttons directly on pencil. This allows the user the option of quickly fine tuning the power output from the field to achieve the desired surgical effect. All pencils are designed with a non-spinning hex collar for multi-positional fail-safe blade locking and a durable, pliable pencil cord to reduce hand fatigue during use and ease connection to the unit.

Also available is a Footswitch Pencil option for those who prefer. Power setting are altered on the base unit, and the pencil is activated by depressing the foot control as required.

Optional Patient Plate is available if required.

Pencil Sheaths

The Sterile Sheath is a convenient and economical way to maintain a sterile field while using the Hyfrecator® electrosurgical reusable pencils. The sheath is designed to fit over both Hyfrecator® and other handheld instruments. The extra sheath length allows the instrument to be placed within the sterile field. A non-sterile sheath is also available.

Employed use of sheaths promotes pencil longevity, speeds changeover from patient to patient, and prevents pencil exposure to contaminants

A non-sterile version is also available to reduce costs where sterility is not required, but protection of the handpiece would be useful.



Electrodes



ConMed Electrolase® disposable electrodes are available in 2 configurations, Sharp for pin-point accuracy and Blunt for general broad-base coagulation. Colour-coding is employed for rapid identification – white for sharp, dark grey for blunt. Disposable Electrolase® electrodes eliminate patient-to-patient cross-contamination and the plastic safety collars simplify handling and help prevent accidental injury.

Available in packs of 50 sterile or 100 unsterile.

Bipolar Accessories

As a multi-function electrosurgical generator, The Hyfrecator® 2000 has a bipolar mode similar to those found in much larger systems.

Bipolar Forceps are designed to grasp, manipulate and coagulate selected tissue. The electric current alternates between the two tips, reducing the length of the current's path and producing a precise, controlled therapeutic effect.

Smooth tipped Adsons pattern Bipolar Forceps are included in the Bipolar



option along with the Foot Switch and Bipolar Cable.

CONMED HYFRECTOR® 2000

CON-001	Hyfrecator®2000 Monopolar & Hand Switching Operation
CON-002	Hyfrecator®2000 Monopolar & Foot Switching Operation
CON-003	Hyfrecator®2000 Monopolar & Bipolar Operation
CON-7-900-1	Telescoping Mobile Hyfrecator® Stand
CON-7-900-5	Hyfrecator® Hand Switching Pencil
CON-7-100-12CS	Electrolase® Sharp Tips - Non-Sterile Box of 100
CON-7-100-8CS	Electrolase® Sharp Tips - Sterile Box of 50
CON-7-101-12CS	Electrolase® Blunt Tips - Non-Sterile Box of 100
CON-7-101-8CS	Electrolase® Blunt Tips - Sterile Box of 50
CON-7-796-18CS	Hyfrecator® Pencil Sheath - Non-Sterile Box of 100
CON-7-796-19CS	Handpiece Sheath - Sterile Box of 25
CON-7-809-11	Hyfrecator® Bipolar Cable
CON-7-809-4	Hyfrecator® Adson Smooth Tips
CON-7-900-4	Hyfrecator® 2000 Foot Switch
CON-7-900-6	Hyfrecator® Foot Switching Pencil
CON-7-900-7	Hyfrecator® 2000 Patient Plate

Other patterns of Bipolar Forceps, plus re-usable Electrodes, are available for special order. Please call for details.

HENRY SCHEIN®
VETERINARY EQUIPMENT



Henry Schein Veterinary Equipment would be delighted to demonstrate this equipment at your practice.

To book a demo please call our Vet Tech Team on 0845 130 9596 or email info@vetinst.com

Animal Stretcher



Can be used by either two people with one at each end or by a single person utilising the wheels at one end. The frame is lightweight coated steel. The Stretcher is made from heavy duty re-inforced nylon. Three 2" (50mm) straps are provided for patient security. 24" (600mm) x 48" (1200mm).

ANIMAL STRETCHER
02742 Animal Stretcher

Soft Carry Stretcher



Without rigid poles the Soft Carry Stretcher can go into tight corners e.g. kennelling and cars. Six padded handles for comfort. Made from heavy duty washable nylon. Easily stored when not in use. Max weight for stretchers is 150lb (67.5kg).

SOFT CARRY STRETCHER
02742A Soft Stretcher 30" (750mm) x 50" (1250mm)
02742B Soft Stretcher 21" (525mm) x 34" (850mm)

Positioning Devices

Redmark Cradles



REDMARK CRADLES
026022 Redmark Cradles - Set of 3 (S, M, L)
026022L Redmark Cradle - Large over 32kg
026022M Redmark Cradle - Medium over 13kg - 32kg
026022S Redmark Cradle - Small up to 13kg

Multi Arm With Limb Brace Attachment In Use



020065



020062




026000

Safety need not be compromised to obtain accurately positioned radiographs. The Multi Arm can replace the human hand in the positioning and holding of limbs for both radiography, surgery and arthroscopy. Also useful for hands free positioning of bronchoscopes & gastroscopes.

- Use cohesive bandage to attach limb
- Three joints allow positioning at any angle.
- A single knob locks all three joints rigid.
- Clamp fixes device to virtually any operating or X-ray table
- Double limb support (for hip radiographs) available as an option.
- As described by Kurt Schulz in Veterinary Surgery

MULTI ARM
026000 Multi Arm - 'Improved' Version includes Single Limb Support
020062 Double Limb Support
020065 Limb Brace Attachment for Multi Arm
MULTIARMSET Multi Arm Set (as above)

Table Top Stand for Hanging Limb Preparation

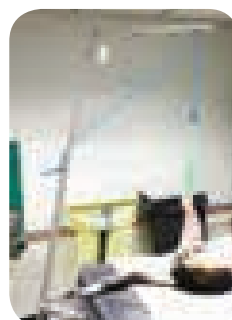

Suspension of a limb can provide significant benefits in both arthroscopy and other orthopaedic procedures. The shoulder joint is distracted and 360 degree access becomes available. Height 1.2m (48") Our limb suspension stand was re-designed to give greater table-top stability. A cleat for tying off limb suspension ropes has been added and the shaft can be dismantled for storage.

TABLETOP STAND FOR HANGING LIMB PREPARATION

R463 Tabletop Stand for Hanging Limb Preparation

Vac Pac Vacuum Positioning Device.



Image courtesy of Pat Ridge

Simply mould the bead filled bag around or under the patient and hook up to any vacuum source, as the air is removed the bag solidifies in the set position. Especially useful for arthroscopy or any surgery requiring secure accurate positioning. The Vac Pac is:

- Re-useable
- Repairable
- Radiolucent
- Warm

VAC PAC
VP23 Vac Pac Size 23 45 x 50cm
VP30 Vac Pac Size 30 72 x 90cm
VP31 Vac Pac Size 31 98 x 90cm
VP32 Vac Pac Size 32 118 x 90cm

Surgical Clippers

Andis Super AGR+ Cordless with Corded Back Up



- * Powerful Clipper designed for all day heavy use.
- * Sensa-Charge Charger recharges battery is less than an hour.
- * High speed motor 3800 strokes per minute.
- * Includes No 10 Ceramic Edge Blade.
- * Cord adaptor converts clipper to mains power in emergencies.

ANDIS SUPER AGR+ CORDLESS

- 0285604 Andis AGR+ Clipper
- 0285142 AGRC Power Groom Cord Adaptor
- 0285605 Andis AGR+ Spare Battery

Oster Clippers A5



The latest version of the industry standard Clipper with a Blade Locking Device to prevent the Blade unclipping in use. (Blade not included). Single speed.

OSTER CLIPPERS A5

- 028080 Oster A5 'Golden' Body Only No Blades

Oster Blades

OSTER BLADES

- 028081 Oster Clipper Blade Size 40
- 028082 Oster Clipper Blade Size 15
- 028083 Oster Clipper Blade Size 50
- 028084 Oster Clipper Blade Size 30
- 028082/10 Oster Clipper Blade Size 10

Cardinal Wet n Dry Clipper Blades

CARDINAL WET N DRY CLIPPER BLADES

- 4412-4 Cardinal Blades x 4

Andis Cordless Pocket Trimmer



Light Duty Trimmer, easily concealed in palm of hand. Handy for i.v. clips or blood tests. Runs off single AA battery (supplied). Cheap enough to have one in every room. Compact ergonomic design.

ANDIS CORDLESS POCKET TRIMMER

- 9861 Andis Cordless Pocket Trimmer Single
- 9861/5 Andis Cordless Pocket Trimmer 5 Pack

Laube 505 Cordless Clipper



A powerful Clipper. Has integrated light for better visualisation.

LAUBE 505 CORDLESS CLIPPER

- 028100 Laube 505 Clipper

Laube Blades

To function with the powerful Laube range of Clippers, Laube have corrected all the faults of all the other blades. Laube Blades fit all Clippers sold by Veterinary Instrumentation, except Pocket Pro Trimmer and Cardinal Wet n Dry.


Highly recommended and very competitively priced.

LAUBE BLADES

- 028111 Laube No. 15 Blade
- 028112 Laube No. 30 Blade
- 028110 Laube No. 40 Blade
- 028113 Laube No. 50 Blade

Nail Clippers

Heavy Duty with Plastic Handles

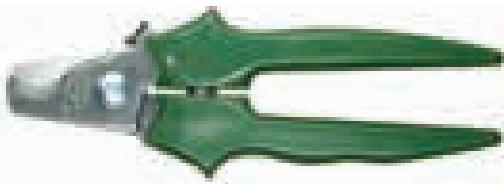

HEAVY DUTY WITH PLASTIC HANDLES

028088 Nail Clipper Heavy Duty Plastic Handle


Compound Action Nail Clippers



COMPOUND ACTION NAIL CLIPPERS

028086 Nail Clipper Compound Action

Nail Clippers Standard

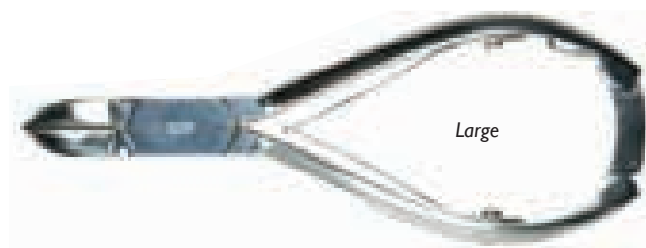


NAIL CLIPPERS STANDARD

028085 Nail Clippers Large 14cm

028085S Nail Clippers Small 12cm

Resco Guillotine Type

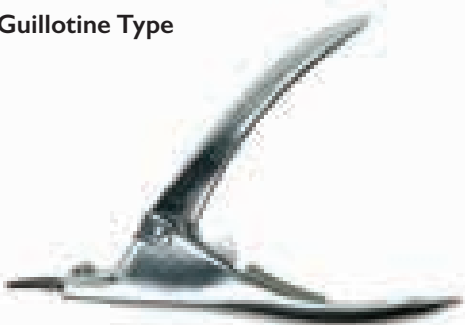

A lightweight clipper suitable for cats and small furries.

RESCO GUILOTINE TYPE

028087 Resco Clippers Guillotine

Beak Trimmer

Based on a familiar mechanism, we have reversed the jaws to allow smooth, controlled and stress free trimming of the lateral extent of the beaks of several species commonly encountered in small animals and exotics practice.

Developed by David Strong



BEAK TRIMMER

028097 Beak Trimmer 65mm

Bird Ring Cutter



Aluminium and plastic rings frequently prevent venous return from the distal limb of small birds. Oedema and necrosis quickly follow unless the ring is removed. Removal of the ring without fracturing the leg is not easy without this special cutting tool.


BIRD RING CUTTER

028090 Bird Ring Cutter

Rabbit Ring Cutters

Identification rings placed on show rabbits can occlude the venous return from the distal hind limb with disastrous results. Removal of the ring without causing additional damage and pain to the limb is very difficult without these purpose designed cutters.

Once the notched jaws are engaged on the ring, firm pressure may be applied to cut the ring without fear of the ring 'flipping over'. Also functions as a general purpose canine nail cutter.

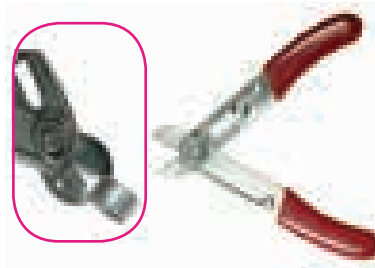

RABBIT RING CUTTERS

028091 Rabbit Ring Cutters

Parrot Ring Cutters



This Device is used by jewellers to remove rings from people's fingers. The ring is cut by turning a hardened saw blade. Parrot rings are very tough and cannot be removed by nail cutters. This Cutter will also remove rabbit rings but bear in mind that the lower blade has to go between skin and ring.

PARROT RING CUTTERS

028093 Parrot Ring Cutter

028094 Spare Blade for Parrot Ring Cutter

Anatomical Models




Canine & Feline Models

When clients better understand their pet's condition they are more likely to accept the treatment you suggest. Models illustrating the common problems are an invaluable aid.

ANATOMICAL MODELS

D1050	Canine Dental Model Clear
D1055	Feline Dental Model Clear
D1057	Leporidae Dental Model Clear
D1062	Set of Canine, Feline & Leporidae Models
AMJ772	Cross-section of Canine Ear
AMJ779	Canine Lumbar Vertebrae with Pathology
AMJ776	Canine Stifle
AMJ775	Canine Hip & Lower Vertebrae
AMJ780	Canine Elbow
AMJ781	Canine Hock

Plastic Training Bones

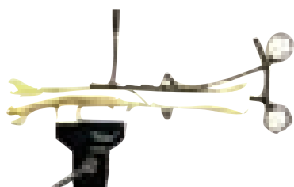

The limited availability of cadaver training programmes means that basic orthopaedic techniques must be practiced on plastic bones. Our Plastic Bones simulate the real thing as closely as possible in that they have cortices and medullary cavities.

PLASTIC TRAINING BONES

2121-20	Femur Right with Block
2117-20	Tibia Right with Block
2117-21	Tibia/ Fibula Right no Block
2003-20	Humerus Right with Block
2002-20	Radius Right Solid Foam with Canal
2002-25	Radius and Ulna Right
2002-22	Radius Solid with Canal/ Transverse Fracture
2003-21	Humerus Foam Cortex with Oblique Fracture
2003-22	Humerus Foam Cortex/ Distal Condylar Fracture
2003-23	Humerus Oblique Condylar Fracture
2003-24	Humerus Growth Plate
2003-25	Humerus High Y Fracture/ Head Fracture
2003-26	Humerus Spiral Fracture
2117	Tibia Left Foam Cortex
2304	Hemi Pelvis Left with Block
2307	Hemi Pelvis Right with Block
2307-1	Hemi Pelvis no Block

Books



A combination of BK02 and BK01 or BK03 will cover nearly all the situations you are likely to face in practice.

BOOKS

BK01	A Guide to Canine and Feline Orthopaedic Surgery by Hamish Denny & Steve Butterworth
BK02	An Atlas of Surgical Approaches the Bones & Joints of the Dog and Cat. 3rd Edition by Donald L. Piermattei
BK03	Handbook of Small Animal Orthopaedics and Fracture Repair. 4th Edition by Donald L. Piermattei & Gretchen L.Flo
BK05	Small Animal Spinal Disorders - Diagnosis and Surgery by Simon J. Wheeler & Nicholas J H Sharp
BK08	Electrocardiography made easy by Larry Tilley
BK09	Feline Orthopaedic Surgery and Musculoskeletal Disease by P M Montavon

DVD/ CD Rom

These DVD/ CD Roms show in detail the following procedures:


- Fracture Fixation including surgical approaches, plate, lag screw and tension band fixation.
- Triple Pelvic Osteotomy using a VI Plate.
- Total Hip Replacement.
- Carpal Arthrodesis.
- Surgical management of Canine CrCL Disease.
- Surgical approaches to the Bones and Joints of the Dog.

DVD/CD ROM

HDFFDVD	DVD Denny Fracture Fixation DVD Version
HDTPDVD	Hamish Denny Triple Pelvic Ostotomy DVD Version
HDTHPDVD	DVD Hamish Denny Total Hip Replacement DVD Version
HDCADVD	DVD Hamish Denny Carpal Arthrodesis DVD Version
DVDSAW	Surgical Approaches DVD Dr Wing Tip Wong
DVDBP	Basic Guide to Bone Plating DVD Dr Wing Tip Wong
DVDALL	Surgical Management of Canine CRCL Disease - covers TTO, TTA, TPLO plus Lateral Suture

Free Training Aids and Brochures


Posters

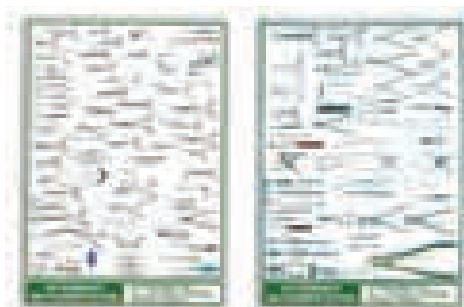

INSTRUMENT POSTERS

- POSTER1** Orthopaedic Poster
POSTER2 General Surgical Instrument Poster

Plate/ TTA Acetates

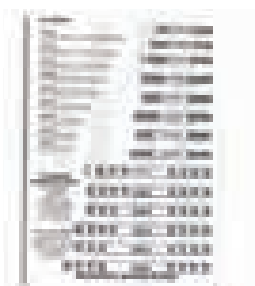


We provide profiles of all our plates on acetate to facilitate correct selection. Simply lay the acetate over the X-Ray. All acetates are to scale – please compensate for X-Ray magnification.

PLATE/ TTA ACETATES

- ACETATES** Acetate Overlay of Bone Plates
TTAPFO Traditional TTA/ MMT Advancement Overlay
TTATAN Common Tangent TTA/ MMT Advancement Overlay


Set by Step Guides

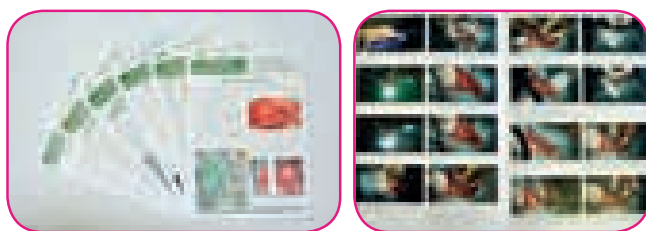

STEP BY STEP GUIDES

- BRTTO** Step by Step Guide to TTO
BRTTA Step by Step Guide to TTA
BRTTAR Step by Step Guide to TTA Rapid
BRTPLO Step by Step Guide to TPLO (Radial Cut)
BRLATSUT Step by Step Guide to Lateral Suture
BRPAT Step by Step Guide to Patellar Luxation
BRLIGAIISO Step by Step Guide to LigaFiba ISO Toggle Sutures
BREXFIX Step by Step Guide to External Fixation
BRTTOVCOT TTO VCOT Paper

Websites

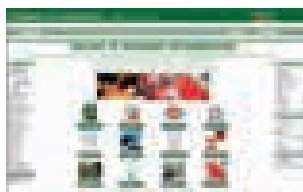

www.vetinst.com

Over 3000 pages of useful information, on-line ordering for over 4000 products, free video downloads and training materials. Over 50 clips of instruments in use.

Catalogues & Brochures

ADDITIONAL CATALOGUES AND BROCHURES


- BR SUTURE** Shoulder Repair Technique Brochure
BR TIGHT Tightrope CCL Technical Brochure
BRAESPOW Aesculap Power System Brochure
BRAESPOWBURBLA Aesculap Power System Burs and Blades Brochure
BRARTHAND Arthrex Instruments Technical Brochure
BRARTHREXEAGLE Arthrex Eagle Eye Camera Brochure
BRCANARTH Canine Arthroscopy Set Brochure
BRCCLS Arthrex Canine Cruciate Ligament Repair Brochure
BRDESOUTLB Desoutter Power Tools Large Bone Brochure
BRDESOUTSM Desoutter Power Tools Small Bone Brochure
BRFIBER Fiberwire Brochure
BRFRVID Dr Fritz Video Compact System
BRLSAS Step By Step Brochure Arthrex Suture Anchor System
BROATS Arthrex Oats System Brochure
BROES Oesophagostomy Procedure Brochure
BRROPES OPES Consol Brochure (electro surgical)
BRSHAV Arthrex Shaver System
BRHELICA Helica Brochure
BRHELISURGMAN Helica Surgical Manual
BRDIXEY Dixey Catalogue

DVDs & CDROMs

DVDS AND CDROM

- DVDALL** DVD Containing: OTT, Lat Suture, TPLO Closing Wedge and Radial Cut, TTA, TTO, Tight Rope and Meniscus Surgery and Full Catalogue
DVDARTHREX Arthrex Multi Brochure DVD
DVDLEHMAN Lehman – Arthroscopy
DVDI/N DVD – Interlocking Nail
DVDBAP Cochran Bisecting Angle Positioner Illustrated CD
DVDTPLOCAD DVD Illustrating Use of Cadmus Plate and Template
CDRAS CD PowerPoint Presentation on Rasenberg Plate (CD ROM)
CDHAU Illustrating Use of Hauptman Retractor

Equine/ Large Animal Catalogue

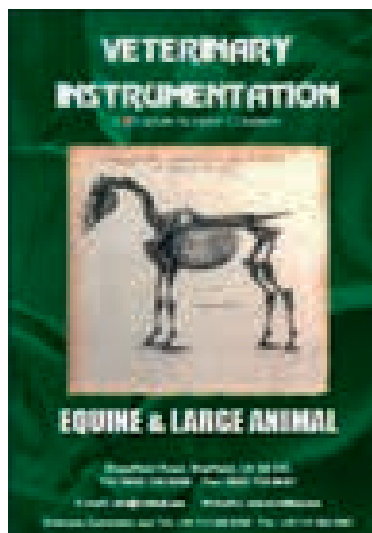

Equine and Large Animal Catalogue with interesting and surgical specialties.

EQUINE/LARGE ANIMAL CATALOGUE

- CATEQ** Equine/ Large Animal Catalogue